



SPECIFICATIONS

Range	T75	0 to 1999 ppm (mg/L)
	T76	0 to 9990 ppm (mg/L)
Resolution	T75	1 ppm (mg/L)
	T76	10 ppm (mg/L)
Accuracy (@20°C)		±2% F.S.
Typical EMC Deviation		±2% F.S.
Temperature Compensation		Automatic, 0 to 50°C with $\beta=2\%/^{\circ}\text{C}$
TDS factor		0.5
Calibration		Manual, at 1 point through trimmer
Probe	T75	MA73075 (replaceable)
	T76	MA73076 (replaceable)
Environment		0 to 50°C; 100% RH
Battery type		3 x 1.5V, alkaline
Battery life		approx. 250 hours of continuous use
Dimensions		165 x 30 x 30 mm
Weight		80 g



ACCESSORIES


MA73075	Replaceable Conductivity probe, LR
MA73076	Replaceable Conductivity probe, HR
M10032B	1382 ppm (mg/L) calibration solution, 20 mL (25 pcs)
M10038B	6.44 ppt (g/L) calibration solution, 20 mL (25 pcs)

ISTR175R1 06/02



USER MANUAL

T75 and T76 Waterproof TDS meters

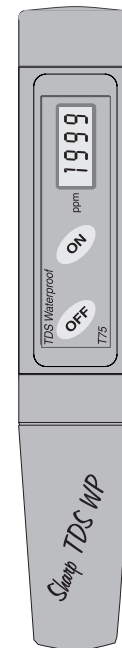

This meter is in compliance with the CE directives.



CALIBRATION:


- Immerse the meter in the proper calibration solution: 1382 ppm (M10032) for T75 and 6.44 ppt (M10038) for T76.
- Unscrew the battery compartment cap on the top of the meter.
- Allow the reading to stabilize and with a small screwdriver turn the calibration trimmer inside the battery compartment until the display shows "1382" for T75 or "6440" for T76.
- Replace the battery compartment cap.



The calibration is now complete.

OPERATION:

- Remove the protective cap.
- Turn the meter on by pressing the ON key.
- Immerse the probe in the solution to be tested.
- Stir gently and wait for the reading to stabilize.
- After use, rinse the electrode with water to minimize contamination.
- Always replace the protective cap after use.

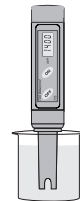


Note: if the reading is over-range, T75 will display "1" while T76 shows "...".



PROBE REPLACEMENT:

The probe can be easily replaced in the following way:

- Remove the protective cap.
- Unscrew the plastic ring on the top of the probe.
- Pull out the MA73075 (or MA73076) probe and replace it with a new one (aligning the colored pin towards the LCD side).
- Make sure the gaskets are in place before screwing back the plastic ring.

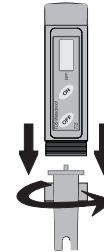

BATTERY REPLACEMENT:

The meter is supplied with BEPS (Battery Error Preventing System) which avoids any erroneous reading due to low battery level, by switching the meter off. To replace the batteries unscrew the battery compartment cap and replace all three 1.5V batteries while paying attention to their polarity. Make sure the gasket is in place before screwing back the cap.

Batteries should only be replaced in a non-hazardous area using the battery type specified in this instruction manual.

WARRANTY:

These instruments are warranted from all defects in materials and manufacturing for a period of **two years** from the date of purchase. The **probes are warranted for a period of one year**.

If during this period, the repair or replacement of parts is required, where the damage is not due to negligence or erroneous operation by the user, please return the parts to either dealer or our office and the repair will be effected free of charge.